

成品尺寸: 120mmX193mm

材 质: 250g铜版纸覆哑膜

RADIANCE XPERT

SIGNIFICANTLY MORE EFFECTIVE THAN SIMILAR PRODUCTS  
ON THE MARKET AT ENHANCING THE APPEARANCE OF FINE LINES

## GOLDEN RADIANCE-8 BEAUTY DEVICE

### 3 MAIN CHANNELS

- Facial rejuvenation
- Targeted skin therapy
- Skin lifting and firming

### 8 MAJOR BEAUTY TECHNOLOGIES

1. Advanced 1-3.3 MHz Multi-Polar Turbo Radio Frequency (RF) Technology: Higher-frequency Turbo RF targets the upper dermis to smooth fine lines and refine texture, while lower frequencies reach deeper layers to restore firmness and boost collagen renewal.
2. 660 nm Red Light Therapy: Promotes skin rejuvenation and enhances radiance for a more youthful appearance.
3. 810 nm Near-Infrared Light Therapy: Enhances absorption of light energy to promote overall skin vitality.
4. 590 nm Yellow Light Therapy: Improves skin tone and texture, reduces redness, and revitalizes dull, tired-looking skin.
5. 415 nm Blue Light Therapy: Calms sensitive, acne-prone skin and helps restore balance for a clearer complexion.
6. 24K Gold Electroplated Treatment Head: Designed to fit facial contours seamlessly, enhancing energy delivery for a more revitalized and refreshed appearance.
7. Multi-Polar Electrical Muscle Stimulation (EMS) Technology: Utilizes impulse and sinusoidal waves to enhance skin firming, toning, and elasticity.
8. Blue LED Charging Base: Utilizes 415nm blue LED light to maintain a clean and hygienic device, ensuring safe and effective treatments every time.

- ◆ The device includes RadianceFirm Advanced Wrinkle Defense & Stem Cell Conductive Gel : Free from parabens, synthetic fragrances, and artificial dyes, this vegan formula is enriched with a blend of powerful peptides, skin-rejuvenating ingredients, and plant-based stem cells to promote youthful, radiant skin.



Use the device **3** times a week.  
The RF and EMS modes each run for **6** minutes, after which the device automatically powers off with a long beep sound.

◆◆◆◆  
For optimal results, dedicate sufficient time to each area of the face during your session. Visible improvements may be seen in as little as 28 days.

## QUICK START GUIDE

### - Product Usage Instructions -

#### 1 Power On/Off and Mode Selection:

- ▶ Short press the power button to turn on the device (you will hear a single "beep" sound).
- ▶ Short press again to cycle through different function modes, with a "beep" sound confirming each switch.
- ▶ Long press the power button to turn off the device (you will hear two "beep" sounds).
- ▶ Blue Light Mode: Balances the skin's appearance, ideal for delicate or acne-prone areas.
- ▶ RF Mode with NIR and Yellow Light: Boosts skin vitality, smooths fine lines, and promotes a radiant, youthful glow.
- ▶ EMS Mode with NIR and Red Light: Enhances skin firmness and delivers a sculpted, toned look through advanced microcurrent technology. Before switching from impulse waveform to sinusoidal waveform, the device emits four "beep" sounds to alert you of the waveform mode change.

#### 2 Application:

- ▶ Wash your face thoroughly prior to using the device. For RF and EMS modes, squeeze the gel onto the treatment head of the device or apply a thin, even layer of gel directly to your face.

**Important:** Ensure that you move the device slowly and continuously over the skin, allowing even distribution of energy across the treatment area. Avoid holding the device in one place to prevent overstimulation and burning.

#### 3 EMS and RF Functions:

- ▶ The RF mode offers two intensity levels: weak and strong, while the EMS mode offers three intensity levels: weak, medium, and strong. Both functions start at the weak mode and can be increased for enhanced results. Adjust the intensity to your comfort level for optimal performance.

#### 4 For optimal results:

- ▶ RF Mode: Gently glide the device in circular motions across each side of your face.
- ▶ EMS Mode: Using gentle upward strokes, glide the device from your jawline toward your cheekbones, to enhance lifting and toning. To target nasolabial folds, start by forming a small 'O' shape with your mouth. Position the device at the tip of your nose, near the corner of your eye. Slowly and gently glide the device down the side of your nose until you reach the corner of your lips. Follow Gua Sha-inspired movements for improved contouring and smoothing.

- 5 After completing one side of your face, repeat the same process on the other side, allowing time for thorough coverage.

- 6 To rejuvenate the skin around your eyes, gently glide the device in soft circular motions while using RF Weak Mode for optimal results.


- 7 After your treatment, leave the gel on your face for 30 minutes before rinsing it off with cool water. Then, apply your preferred skincare products. To clean the device, gently wipe the RF treatment head with a damp cloth, wet wipes, or paper towel.




Forehead lines



Under-eye line / Crow's feet

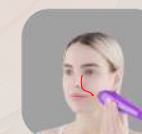

Nasolabial folds

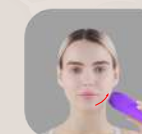

Jawline

**Note:** Always remember to avoid soaking the device in water.

**Disclaimer:** It is your responsibility to use the device properly. Follow the instructions carefully to avoid burns or injury. This product is intended for cosmetic use only.